


APSW 4-PLY, WIRE-REINFORCED, MANDREL-WRAPPED, SILICONE SUCTION HOSE

PLATINUM-CURED; SMOOTH OR CONVOLUTED O.D.

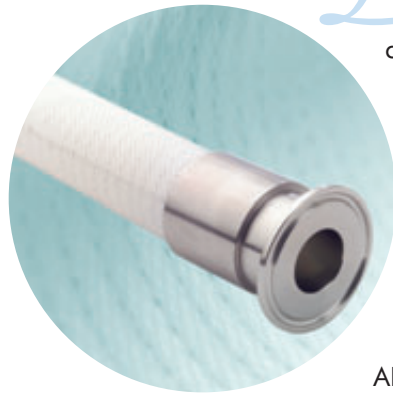

Low volatile grade, platinum-cured silicone hose, produced in our certified **Class 1000 clean room** for critical pharmaceutical, biomedical, cosmetic, and food applications, comprises the core of APSW. The core is then slipped onto a mandrel and wrapped with polyester mesh fabric, 316 stainless steel reinforcing wire, and additional silicone. It is cured into a homogeneous hose to enhance its pressure and vacuum capabilities. APSW has undergone extensive physical, chemical, and biological testing and meets USP **Class VI**, FDA CFR 177.2600, ISO 10993, European Pharmacopoeia 3.1.9, and 3-A standards.

APSW is not intended for implantation and is not to be used for continuous steam applications.



SILICONE

www.advantapure.com

KEY FEATURES



Core

- Available with either a smooth or convoluted O.D.
- Low volatile grade silicone suitable for pharmaceutical, biomedical, cosmetic, and food applications such as ultrapure water transfer, food processing, washdown stations, skid systems, and other suction and discharge uses
- Core made of **NSF 51** listed silicone (National Sanitation Foundation)
- Hardness value of 70 Shore A
- Rated for full vacuum
- Crush and kink resistant
- Durable and heavy-duty yet translucent for visual flow contact
- Two styles available: Smooth O.D. or, for an improved range of flexibility, Convoluted
- General temperature range: -100°F (-73.3°C) to 400°F (204.4°C)
- Sterilizable by autoclave, CIP, SIP, and gamma radiation processes
- Resists temperature extremes, chemical attack, ozone, radiation, moisture, and environmental exposure
- Imparts no taste or odor to critical streams
- Stocked in 12 ft. straight lengths
- Manufactured under strict Good Manufacturing Procedures in a tightly-controlled environment
- Wrapped in clear polyethylene at manufacture to ensure cleanliness
- Documented lot traceable with identification on bags and boxes
- Documented quality control
- Complete validation package available upon request
- Custom lengths, sizes, colors, special cleaning, and/or packaging available (lead times and minimum order quantities may apply – call for details)
- Available with Nomex® braid construction – call for details

