

# MOLDED SILICONE ASSEMBLIES

## PLATINUM-CURED SILICONE TUBING AND FITTINGS



Supplied as silicone fittings attached directly to silicone tubing or hose, AdvantaPure's Molded Silicone Assemblies are well-suited for applications involving aseptic transfer, sterile processing, and other critical applications in the pharmaceutical, biomedical, and high purity industries. The assemblies reduce installation time, improve cleanliness, and retain the many benefits of silicone.

Molded Silicone Assemblies are not intended for implantation or continuous steam applications.



### KEY FEATURES



- Eliminates entrapment issues and leaks associated with barb-type fittings
- Provides a smooth inner surface with a seamless transition
- Offers consistent interior dimensions for maximum flow rates
- Reduces end-user assembly times—arrives ready to install
- No product contact with metallic materials
- Retains silicone's temperature resistance: -100°F to 400°F
- Extremely flexible
- Sterilizable by CIP, SIP, and autoclave processes
- Imparts no taste or odor to critical streams
- White-room manufactured
- Sealed in clear polybags during manufacturing process to ensure cleanliness
- Documented lot traceable with identification on bags and boxes
- Documented quality control

SILICONE

[www.advantapure.com](http://www.advantapure.com)


# MOLDED SILICONE ASSEMBLIES

## PLATINUM-CURED SILICONE TUBING AND FITTINGS

Molded Silicone Assemblies are available from AdvantaPure in mini and standard Tri-Clamp® fitting styles and are supplied with male gaskets molded directly to the face of the fitting. Straight-through and step-up bore configurations are produced in sizes up to 1" tube I.D. Protective backup cups, offered in polysulfone or series 304 stainless steel

for extra durability, provide a stable clamping surface and safeguard the fittings during installation and use.

All assemblies are custom made to best meet each application's requirements and incorporate either unreinforced silicone tubing or braid-reinforced silicone hose. Working pressure ratings vary depending on product. Call AdvantaPure today to discuss your needs.

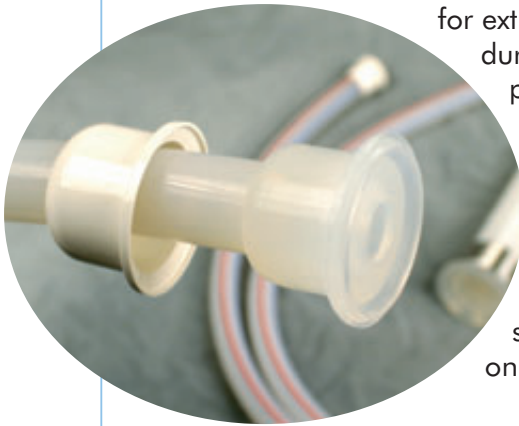

SILICONE

## Purity in Fluid Flow Systems<sup>SM</sup>

Molded Silicone Assemblies are part of a complete line of products for the Pharmaceutical, Biomedical, Cosmetic, and Food industries. Items available include . . .



- Unreinforced, 50 Shore A, platinum-cured silicone tubing
- Braid-reinforced, platinum-cured silicone hose
- Wire and mesh reinforced, mandrel-wrapped, platinum-cured silicone suction hose
- Stainless steel overbraided PTFE hose; rubber-covered FEP and EPDM hose
- Sanitary fittings and complete hose assemblies including molded silicone ends
- Identification solutions including the Hose Track™ Process Equipment ID & Lifecycle Analysis System
- Silicone injection molding capabilities – manifolds, stoppers, and sealing systems for single use processes

Address: NewAge® Industries AdvantaPure®  
145 James Way  
Southampton, PA 18966 U.S.A.

**AdvantaPure**  
NEWAGE INDUSTRIES



DISTRIBUTED BY:

Toll Free Phone: **888-755-4370**  
Phone: **215-526-2151**  
Toll Free Fax: **888-258-4293**  
Fax: **215-526-2167**  
Web Site: **www.advantapure.com**  
E-mail: **sales@advantapure.com**

**ProcessHQ, Inc.**  
21935 Plummer Street  
Chatsworth, CA 91311  
818 341-1782 (voice)  
818 475-1589 (fax)  
info@processHQ.com

©2007 NewAge® Industries, Inc.  
AdvantaPure products are produced by NewAge® Industries, Inc.  
Hose Track™, NewAge® Industries AdvantaPure®, the AdvantaPure logo, NewAge® and the "N" logo are registered trademarks of NewAge® Industries, Inc. • Purity in Fluid Flow Systems™ is a service mark of NewAge® Industries, Inc. • Tri-Clamp® registered trademark of the Ladish Co.



**ESOP** Employee Owned  
for Your Benefit

Celebrating over  
**50**  
Years  
of Service  
1954-2007

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term NewAge.